




When your home is your castle, the thought of moving can be heart breaking.

But if it's the stairs that are presenting a problem, then a stairlift may be the answer to **regaining your independence**, and making your home fully accessible again.

National Stairlifts can customise a safe, and comfortable stairlift solution that is affordable and easy for you to use.

We are a family-owned business with years of experience in all types of stairlifts. Our team can custom design a solution to suit your physical requirements, your stairs, and your budget.

Our stairlifts meet Australian Standard 1735.7 and are design registered in Australia.

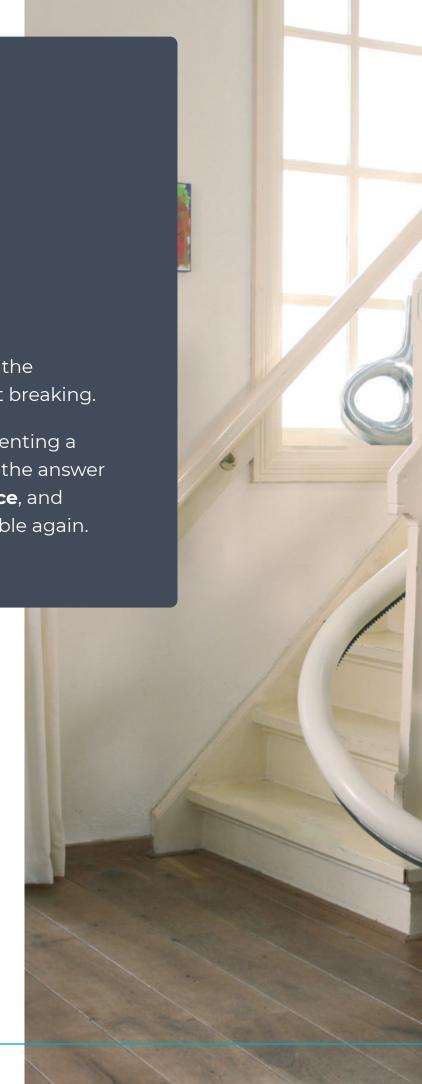



The freedom to live with independence in your own home.

National Stairlifts is proud to be a partner of Handicare, a market leader in the field of home care, with a reputation for the manufacture of innovative, quality stairlifts. Handicare stairlifts are favoured by people around the world, because they are safe, affordable, and easy-to-use. With a choice of colours and finishes, the stairlifts also integrate seamlessly with many home décor schemes.

Our range of stairlifts suit both straight and curved staircases, for both indoor and outdoor applications. Each track is customised to suit its setting, and the chairs are tailored to suit the user's needs.



WIDTH UNFOLDED

HANDICARE 1000

INDOOR USE

The Handicare 1000 is one of the slimmest straight stairlift tracks on the market. With no gearing exposed, the 1000 track offers a sleek look and takes up little room along your staircase to ensure others can safely get by.

- Fold-up seat, footrest and armrests save space for other stair users
- Safety sensors immediately stop the lift if anything is detected in its path
- Adjustable seat height ensures the perfect fit for each user
- Manual swivel seat for safer transfer onto and off the lift (powered upgrade option available)

- Key switch prevents unauthorised use
- Vinyl seat for easy cleaning
- Anodized aluminum track with no exposed gearing
- Large, carpeted footrest for greater traction
- Manual emergency stop
- Two remote controls included
- Battery operated the lift will still run in the event of a power failure





MAXIMUM LOAD 200kg



CHARGE SUPPLY
Multiple charge
points



WIDTH UNFOLDED 630mm



WIDTH FOLDED 350mm



MIN STAIR WIDTH 655mm



SEAT COLOUR Cocoa, ruby, sand, sapphire, jade, slate

HANDICARE 1000

OUTDOOR USE

The Handicare 1000 is one of the slimmest straight stairlift tracks on the market. The outdoor stairlift is UV protected and includes a waterproof cover so it is always ready for use.

- Fold-up seat, footrest and armrests save space for other stair users
- Safety sensors immediately stop the lift if anything is detected in its path
- Adjustable seat height ensures the perfect fit for each user
- Manual swivel seat for safer transfer onto and off the lift

- Key switch prevents unauthorised use
- Vinyl seat for easy cleaning
- Anodised aluminum track with no exposed gearing
- · Manual emergency stop
- Two remote controls included





MAXIMUM LOAD 160kg



CHARGE SUPPLY
Multiple charge
points



WIDTH UNFOLDED 630mm



WIDTH FOLDED 350mm



MIN STAIR WIDTH



SEAT COLOUR

HANDICARE FREECURVE

INDOOR USE

The Freecurve is a revolutionary design which utilises a single, custom tube built to fit your staircase. Coupled with several proprietary design options, the Freecurve maximises your comfort while providing optimal aesthetics and convenience.

- LED illuminated footrest improves visibility in the dark, reducing the risk of stair users from stumbling over the unfolded footrest
- Includes a user-friendly retractable seatbelt
- Includes two handheld remotes to easily send and retrieve the stairlift
- Key operation prevents unauthorised use
- Includes a variety of seat options and track colours, and the ability to custom colour match your track to your interior





MAXIMUM LOAD 125kg



CHARGE SUPPLY Multiple charge points



WIDTH UNFOLDED 640mm



WIDTH FOLDED 420mm



MIN STAIR WIDTH



SEAT COLOUR Beige, anthracite, red, cream, cocoa brown

HANDICARE SIMPLICITY 2000

OUTDOOR USE

Made with weather-resistant materials, the Simplicity 2000 Outdoor comes with a waterproof cover to protect it from the elements so it is always ready for use.

- Weight limit: 115kg or 137kgs
 with heavy duty kit
- Fold-up seat, footrest and armrests save space for other stair users
- Safety sensors immediately stop the lift if anything is detected in its path
- Adjustable seat height ensures the perfect fit for each user
- Manual swivel seat for safer transfer onto and off the lift

- Soft start/stop
- Key switch prevents unauthorised use
- · Vinyl seat for easy cleaning
- · Manual emergency stop
- Two remote controls
- · Retractable seatbelt
- Battery backup in the event of power outage





MAXIMUM LOAD 137kg



CHARGE SUPPLY
Multiple charge
points



WIDTH UNFOLDED 690mm



WIDTH FOLDED 400mm



MIN STAIR WIDTH 720mm



SEAT COLOUR Grey

Stairlift Options

Our stairlifts afford a variety of options with different levels of automation and features to meet your physical needs, personal preferences and budget. Choose from a range of options to create a stairlift that is right for you and your home.



Restraints

We offer a range of seat belts (clip, velcro and retractable), body harnesses, and foot restraints.




Powered Swivel Seat

When you arrive at the top of the stairs you will need to swivel your seat around to exit onto the landing safely. If you choose the powered swivel option, the seat will automatically swivel when you reach the top of the stairs.



Powered Footplate

If you choose the powered footplate option, it will remove the need for you to bend down to fold or unfold the footplate.

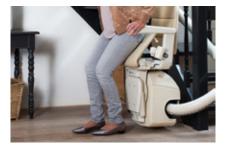


Active Seat

The active seat provides assistance for those users who would like some additional support when getting on or off the stairlift. It can also provide a greater sense of security when sitting down as the seat can be felt directly behind you.



Additional Seat Design

Users who have difficulty bending at the hip or knee joint may choose a 'Perch Seat' – this option reduces the angle of the sitting position and therefore can reduce some people's discomfort. We also provide options for children's seats.

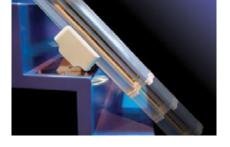

Slide Track

If there is a walkway at the top or bottom of the staircase, a slide track is useful to move the rail and the stairlift out of the way.



Hinged Track

If there is a walkway at the bottom of the stairs, a hinged track will fold out of the way.

Visit our showroom today, call us on **1800 940 575** or visit **nationalstairlifts.com.au** for more information.



How do I choose the right stairlift?

To ensure that your stairlift suits your needs, one of our team can visit your home to assess your staircase and discuss your requirements. This short consultation is without any obligation, and will be conducted by one of our qualified installers.



Each stairlift is designed to suit your individual staircase.

Using our specialist PhotoSurvey software we take measurements of your staircase, and then create an accurate 3D model – this software ensures the stairlift track is designed to perfectly fit your space. We also advise of any adjustments required to fit your stairlift, for example moving a heater, and which side of the staircase would best suit your stairlift.



What are my stairlift options?

Stairlifts are either manufactured for straight or curved staircases, and for internal and external applications. Each of these stairlifts provide a range of features and functionality, and depending on your needs and budget there are options to increase the automation of the lift to make it even easier to use.

Visit our showroom today, call us on **1800 940 575** or visit **nationalstairlifts.com.au** for more information.

This publication provides general information only. No representation or warranties are made regarding the accuracy, completeness or relevance of the information provided. All measurements quoted are approximate values only and can vary. National Stairlifts accepts no liability for any loss or damage resulting from any reliance on the information provided in this publication. National Stairlifts is a registered business name.

© January 2020 National Stairlifts Pty Ltd ABN 88627357386. All rights reserved. This guide or any part of it may not be reproduced without prior written consent from National Stairlifts Pty Ltd. 'I support NDIS' and 'Registered NDIS Provider' trademarks owned by the Commonwealth Government.

QBCC 15102469

QLD Electrical Licence 85362 NSW Electrical Licence 331996C VIC Electrical Licence 28727 ACT Electrical Licence 20181244 SA Electrical Licence 293661 NT Electrical Licence C3911



GET IN TOUCH

T 1800 940 575
E info@nationalstairlifts.com.au
W nationalstairlifts.com.au

ABN 88627357386

National Stairlifts Pty Ltd

15/40 Counihan Rd Rocks Riverside Industry Park Seventeen Mile Rocks QLD 4073